

Independence with Ergo Comfort

Ergo Traveller

KP-10.3S



The new Karma “Ergo Traveller” is a lightweight, compact, manoeuvrable power chair which now incorporates Karma’s patented ‘S-Ergo’ seating system. It is compact enough for indoor use but powerful enough for outdoor use. The concept behind the design of the power chair was to make it as light as possible and convenient to transport.

The “Ergo Traveller” can be taken apart quickly and easily into four manageable components. The heaviest part weighs just 24kg. The quick release battery box has a convenient charging point built into it, which allows the batteries to be charged away from the power chair.

Why not try one at your local Karma retailer?

S-Ergo



Pressure redistribution





Reduced sliding

Traditional seat



Increased risk of sores due to load on pressure points.

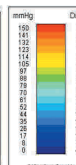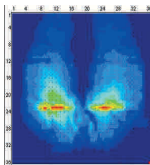

S-Ergo seat



The S-Ergo system (intelligent s-shaped ergonomic seating) provides efficient pressure relief by spreading weight over a greater area, at the same time provides stability and reduces sliding.



Can be taken apart quickly and easily.

Model															
KP-10.3	103	62	94	43 17"	F: 200x50mm R: 12.5"x2.25"	200Wx2	22Ahx2	VR 2	82.5	8°	6km	16-20km	Steel	115kg	49kg

▶ The above are estimated data under 20°C (68°F), an 80kg (176lbs) user, brand new fully charged batteries and a constant driving speed of 6km/hr (3.7mph). Distance, climbing ability and operating time can vary according to condition. ±1cm ±0.5kg

※ KARMA reserves the right to modify information herein without further notice.

Features

- Ergonomic seat system provides comfortable seating.
- Can be taken apart in seconds.
- Quick-release battery pack.
- Wipeable upholstery.
- Height adjustable and flip back armrest.
- Adjustable backrest angle.
- height adjustable footplates.
- Right or left hand joystick control positions
- Heaviest component just 24kg.



Quick release battery



Compact to fit your boot



Height adjustable armrest



Adjustable back angle



Height adjustable footrest

Sold and Serviced by:

104-40653 Rev.00 SEP.2011

• KARMA reserves the right to modify information herein without further notice. •